

**Technimount System** 

Administration
445 Papin Street—Office 103
Quebec City (Quebec) Canada
G1P 3T8

Telephone: +1-581-998-9820 Fax: 1-855-339-6351

info@technimount.com www.technimount.com

## Bracket Pro Serie™ 25-WB

Mounting Bracket system for Monitor/Defibrillator ZOLL® X Series®

Any medical equipment that has been inadequately installed in an ambulance, other emergency medical service (EMS) vehicle or crash cart, can potentially pose a serious risk or become dangerous in

### **Features**

- > SAEJ3043 certified
- Lightweight
- Compact size
- Designed for monitor without bag
- Swivel base 360° with locking system
- Quick Release and reattachment
- Safe and Easy
- Immediate access
- Match with other mounting systems by Technimount

#### **Application of use**


- Mainly for hospital environments
- Crash Cart
- Ambulance
- EMS Vehicle
- > Stretcher/Cot
- Military Vehicle

#### **Advantages**

- Minimum space for installation
- No vertical release
- Biological fluid and sodium resistant
- Rust resistant
- 1-year full warranty

the event that the vehicle makes frequent or sudden stops, turns or accelerates quickly. This situation can result in damage to expensive medical equipment or even inflict injury on a patient, ambulance personnel or clinic staff.

The \*Bracket Pro Serie\*\*\* 25–WB (Without bag)\* is designed to give fast and easy access to the unit during emergency use. Features include small footprint design, lightweigh

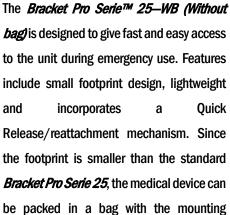

system. Being lightweight, the end user won't even notice the mounting bracket when attached to a medical device.

To spare the expenses of costly repair, replacement of the equipment, and possible legal redress for injury, consider outfitting your EMS vehicle or crash cart with a deployment system for monitors and/or defibrillators. Our mounting systems offer different configurations for most external defibrillators on the market. The deployment of the unit is almost instantaneously when engaging the release button on the base. The complete assembly has been tested for "G" Force tolerances and compliant with SAE standards.

#### **About Us**

**Technimount System** creates mounts and brackets for portable medical devices such as defibrillators, ventilators, and pumps, that adapt to different brands of medical equipment for OEM ambulance manufacturers, EMS providers, and hospitals.

For use with model **ZOLL®** Defibrillator "X Series®"

# **TECHNICAL SPECIFICATIONS**

Product name:	Bracket Pro Serie™ 25—WB
Compatibility model:	Monitor—Defibrillator ZOLL® "X Series®" or "Propaq MD"
Working Environment:	EMS/Hospital—Clinic—BioMed Dept. —Military hospital
Material:	Base: Aluminum— high density w/anodize process
	Device section: Aluminum— high density w/anodized process
Temperature of use:	-35° C@ 45° C (Direct conduction)
Impact Resistance1:	(30 "G") Plus 10% safety factor SAE J3043 CRIQ ref.48064
Installation points:	4 points on the surface base: 4 points on the device mount
Dimensions (Assy—overall):	13-1/4" Width X 9-3/8" Depth <b>X</b> 2-1/5" Thickness (inch).
(Bracket Only):	6-1/2" Width X 6-5/16" Depth X 1-1/4"
Rotation:	360 degrees with locking system (Standard Surface Base)
	360 degrees with locking system (Extended Surface Base)
Features:	With or without carry case protector
	*Compatible with <b>Technimount Stretcher Safety Arm System</b>
	*Compatible with Extended Surface Base
	*Free access to all connectors and accessories (top)
Weight:	Total (2 parts) 2.9 lbs.—device section: 1.35 lbs.
Options:	Extended Surface Base (longer)
	Spacer Pro (height extension mount 3/8 ", ½", ¾", 1",1 ¼")
	Converter Pro (surface base cover for work space)
Locking mechanism resistance:	Tested 20,000+ "press & release" activations w/o any reduction in use
Part number <sup>2</sup> :	100-00-XZ-WB "Bracket Pro Serie™ 25—WB" Complete assembly
Part number	100-11-XZ-WB "Bracket Pro Serie™ 25—WB" Bracket Device Only
Part number	100-20-UN "Bracket Pro Serie™ 25" Surface Standard Base Part Only





<sup>&</sup>lt;sup>1</sup> *Technimount System* works to meet and exceed all legal requirements in the deployment system for EMS vehicles, according to the Provincial & Federal Government regulations. It may not meet your local requirements, please contact your administrator.

<sup>&</sup>lt;sup>2</sup> Part number must be confirmed prior sending purchase order. Technimount System reserves the right to make changes on the product, documents and technical specifications without notice and responsibilities. DOCUMENTS PRINTED IN CANADA CGC Technologies – Division of medical equipment & accessories – Registered brand 2013 ZOLL® and X Series® are trademarks and registered trademarks of ZOLL Medical