




**TECHNIMOUNT**  
EMS™

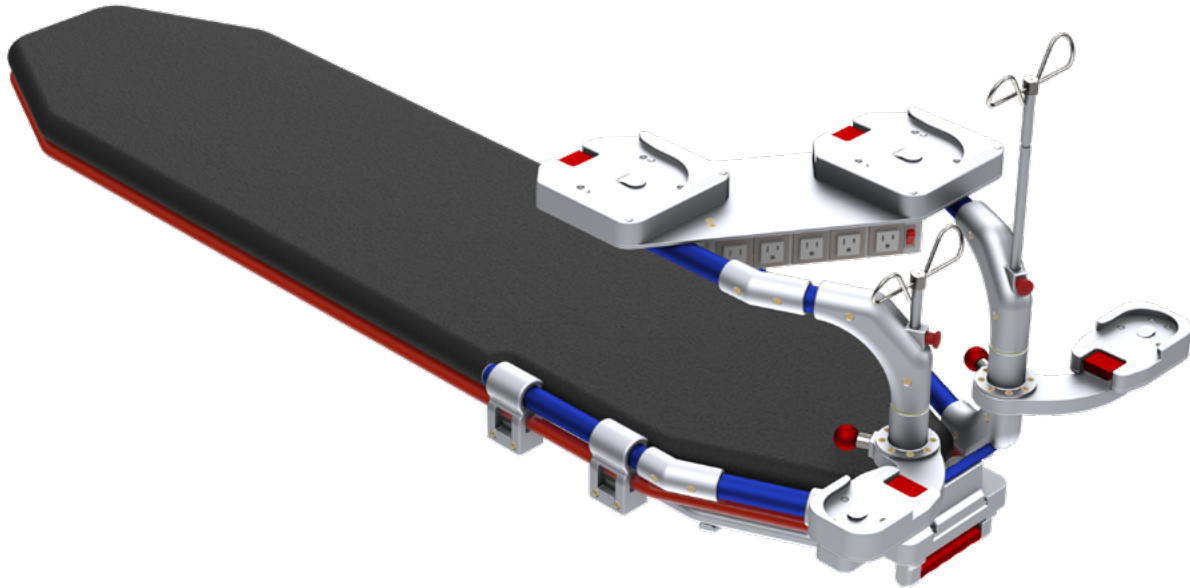

## **XTENSION PRO™ ASSISTANT - LP**

MOUNTING SYSTEM FOR AIR TRANSPORT OF  
MULTIPLE MEDICAL DEVICES ON THE LIFEPORT™ COT

### **SAFETY & FLEXIBILITY WHERE IT MATTERS MOST**

The Xtension Pro™ Assistant - LP is designed for air critical care transport. It facilitates lateral patient transfer from bed to cot and improves accessibility to patients during transport.

The Xtension Pro™ Assistant - LP is tested in compliance with FAA, 14 CFR § 23, ensuring the safety of patients and crews during air transport.



**Air**

At Technimount EMS, we understand that EMS and Critical Care Transport teams need solutions that they can rely on. We design and develop quality mounting systems which are safe and flexible, allowing you to concentrate on what matters most - patient care.


# XTENSION PRO™ ASSISTANT - LP

## POTENTIAL CONFIGURATIONS


Monitor/Defibrillator



Ventilator

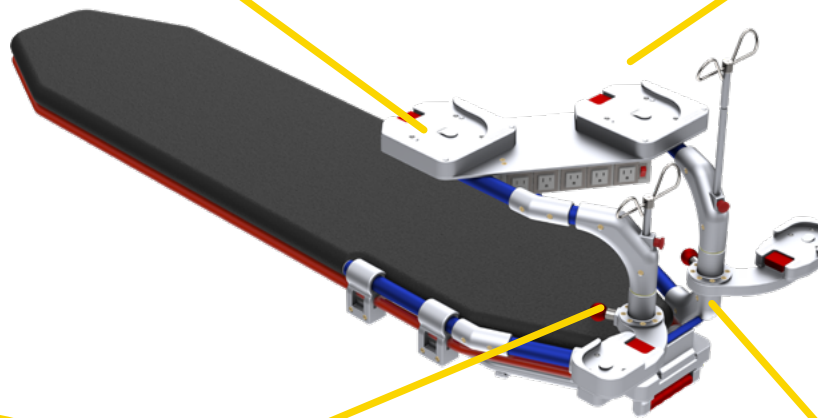

Pumps

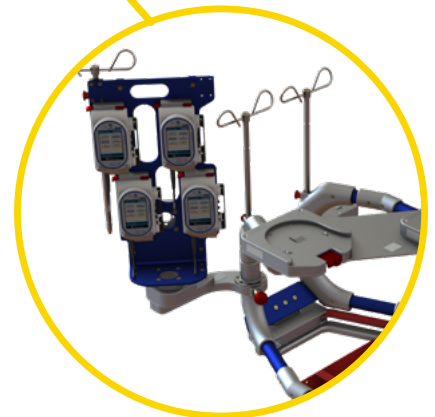


Pumps

# TECHNICAL SPECIFICATIONS

<b>Product Name</b>	Xtension Pro™ Assistant - LP
<b>Product Description</b>	Mounting System for Air Critical Care Transport of Multiple Medical Devices on the LifePort Cot.
<b>Compatibility of Stretcher*</b>	LifePort Cot
<b>Compatibility of Mounting Systems</b>	Compatible with Technimount's: <ul style="list-style-type: none"><li>• Bracket Pro Serie® Mounting Systems with Standard Bottom Disc and Bottom Micro Disc</li></ul>
<b>Environment of Usage</b>	CCT (Flight)
<b>Product Section Material</b>	<ul style="list-style-type: none"><li>• Main Assembly and Clamp Blocks and Locking Handle: Aluminum 6061-T6</li><li>• Cot Support Structure: Aluminum 5052-H32</li></ul>
<b>Product Finish</b>	<ul style="list-style-type: none"><li>• Main Assembly: Blue-Anodized Finish</li><li>• Cot Support Structure, Equipment Support Structure and Clamp Blocks: Silver-Anodized Finish</li></ul>
<b>Environmental Regulations</b>	Complying with RoHs 3
<b>Temperature of Usage</b>	31°F to 113° F (-35° C to 45°C)
<b>Certification</b>	Tested in Compliance with: <ul style="list-style-type: none"><li>• FAA, 14 CFR § 23 (applicable sections)</li></ul>
<b>Installation</b>	<ul style="list-style-type: none"><li>• Locking System - Main Assembly inserts into 4 Clamp Blocks Installed on the Cot</li><li>• Locking Handle located at the Foot End of the Cot is used to Easily Remove the System from the Cot</li></ul>
<b>Features</b>	<ul style="list-style-type: none"><li>• 2 Standard Surface Bases with 360° Rotation</li><li>• 2 Swivel Micro Shelf with 4 Different Positions</li><li>• Mobile Equipment Support Structure</li><li>• 2 integrated IV Poles</li></ul>
<b>Positioning Configuration</b>	Optional mirror positioning for the upper tray - Available upon request
<b>System Physical Overall Dimensions</b>	Width 19.0 in (48.3 cm), Depth 28.835 (73.3 cm), Height 20.577 (52.3 cm)
<b>Mobile Equipment Support Structure</b>	Width 19.0 in (48.3 cm), Depth 25.4 in (72.1 cm), Height 18.6 in (47.3 cm)
<b>Mobile Equipment Support Structure Complete System Weight (With brackets, Without Medical Devices)</b>	<ul style="list-style-type: none"><li>• 52.43 lb (132.08 kg)</li><li>• 68.02 lb (172.72 kg)</li></ul>
<b>Maximum Estimated Loading Capacity</b>	66.0 lb (29.94 kg)
<b>Part Number</b>	1610-00-LPT01
<b>Other Models and Configuration Options</b>	Contact Technimount EMS for more information or for other options at <a href="mailto:customerservice@technimount.com">customerservice@technimount.com</a>

\* Product and medical equipment manufacturers names are trademarks™ or registered trademarks® of their respective holders. Technimount does not have a commercial relationship with these medical equipment manufacturers.

# YOUR SAFETY IS OUR PRIORITY

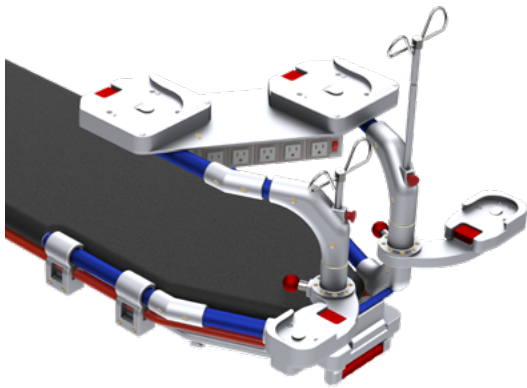

## FEATURES

- ✓ Secures multiple medical devices during Air Critical Care Transport
- ✓ Two Standard Surface Bases (360° rotation) compatible with Technimount's Bracket Pro Serie with Standard Bottom Disc
- ✓ Two Swivel Micro Shelf with four different positions compatible with Technimount's Bracket Pro Serie with Bottom Micro Disc
- ✓ Locking handle located at the foot end of the cot
- ✓ Made from high-density aluminum with blue and silver-anodized finish
- ✓ Tested in compliance with FAA, 14 CFR § 23 (applicable sections)

## BENEFITS

- + Open-sided and tailored design facilitates lateral patient transfer from bed to cot and optimizes space within the aircraft
- + Improves medical equipment management during critical care transport
- + Increases flexibility and access to medical devices from all angles
- + Improves accessibility to patient during air transport
- + Certified mounting system that complies with the highest industry standards and regulations
- + Reduces the risk of damage to the equipment and ensures the safety of patients and crews during air transport



## SAFETY AND FLEXIBILITY WHERE IT MATTERS MOST

[technimount.com](https://www.technimount.com)

