




**TECHNIMOUNT**  
EMS™

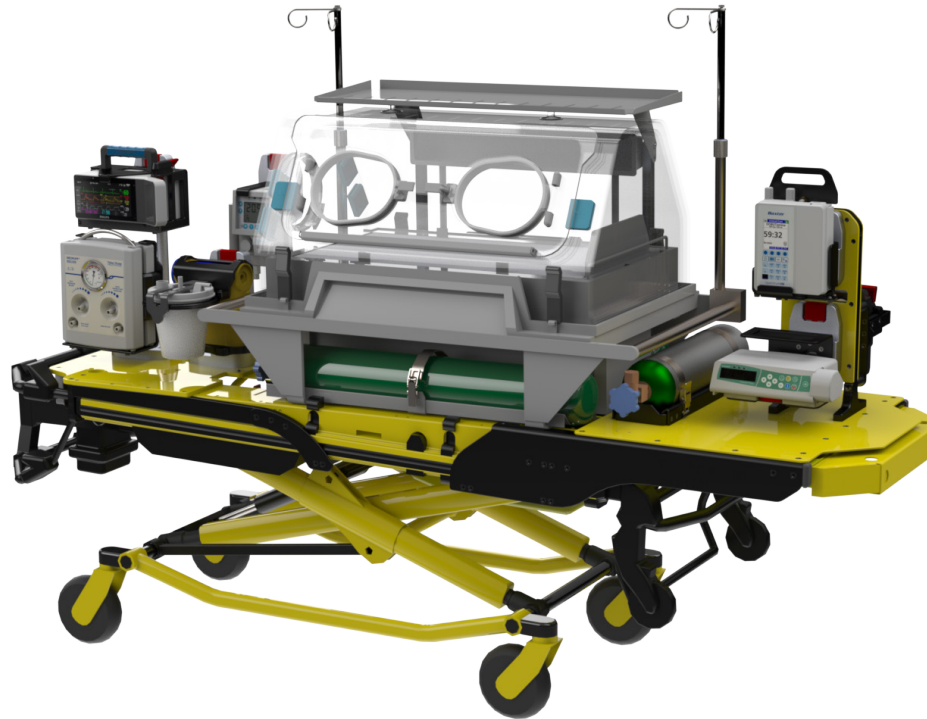

## NEONATAL STRETCHER SYSTEM

MOUNTING SYSTEM FOR THE STRYKER POWER-PRO™ IT STRETCHER

### SAFETY AND FLEXIBILITY WHERE IT MATTERS THE MOST

Our solution is designed for the safe transport of medical devices during neonatal inter-hospital transfers in critical care settings.

This scalable system ensure the centralization of the necessary equipment on the stretcher by allowing the optimization of patient transport.



Ambulance



Stretcher



Hospital

At Technimount EMS, we understand that EMS and Critical Care Transport teams need solutions that they can rely on. We design and develop quality mounting systems which are safe and flexible, allowing them to concentrate on what matters most, patient care.

# NEONATAL STRETCHER SYSTEM

## CONFIGURATION OF THE SUPPORT STRUCTURE

### STRUCTURE COMPONENTS

Patient Foot  
End Support  
Structure

Patient Head  
Support  
Structure




### ASSEMBLED STRUCTURE ON THE STRETCHER



# NEONATAL STRETCHER SYSTEM

## AVAILABLE SETTINGS

Safety MD-Transporter  
Standard and Short Version  
(allows rotation of medical device)

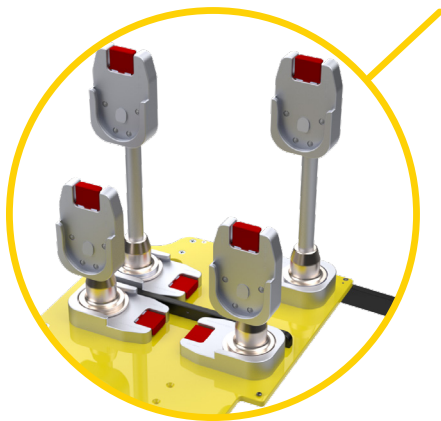

Medical grade power bar with  
surge protection (6 outlets)




Battery Charger Bracket



Vertical Mounting System



Gas Cylinder Bracket

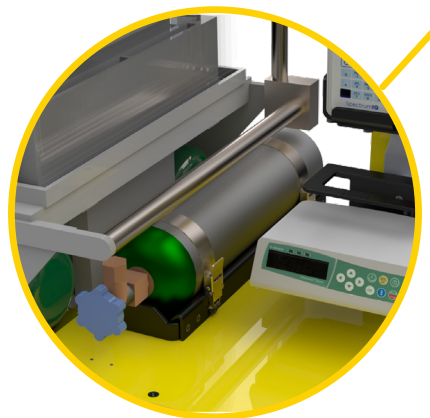

Support Structure





# NEONATAL STRETCHER SYSTEM

## POTENTIAL CONFIGURATIONS

Maxtec MaxBlend Low Flow™  
Oxygen Blender Bracket

Philips IntelliVue™ X3  
Monitor Bracket

Baxter Spectrum IQ  
Infusion System Bracket

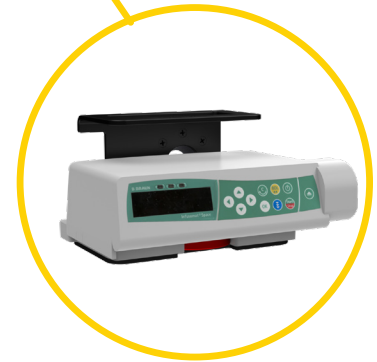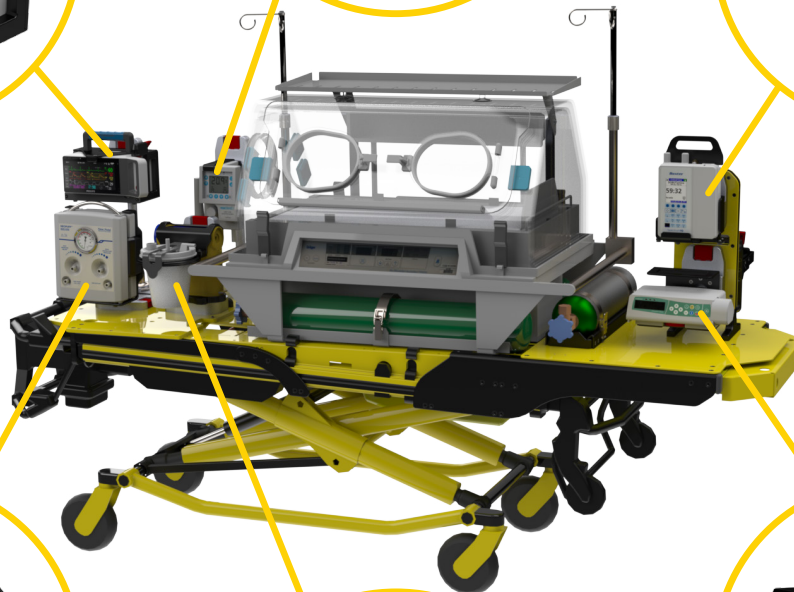

NEOPUFF™ Infant T-Piece  
Resuscitator Bracket

Laerdal Compact  
Suction Unit® 4 Bracket

B. Braun Infusomat® Space/  
Perfusor® Space Pump Bracket

# NEONATAL STRETCHER SYSTEM

## POTENTIAL CONFIGURATIONS


Maxtec MaxBlend™ 2  
Oxygen Blender Bracket



ICU Medical ESP 60™  
Syringe Pump Bracket



Bio-Med Devices 2003FL™  
Oxygen Blender Bracket

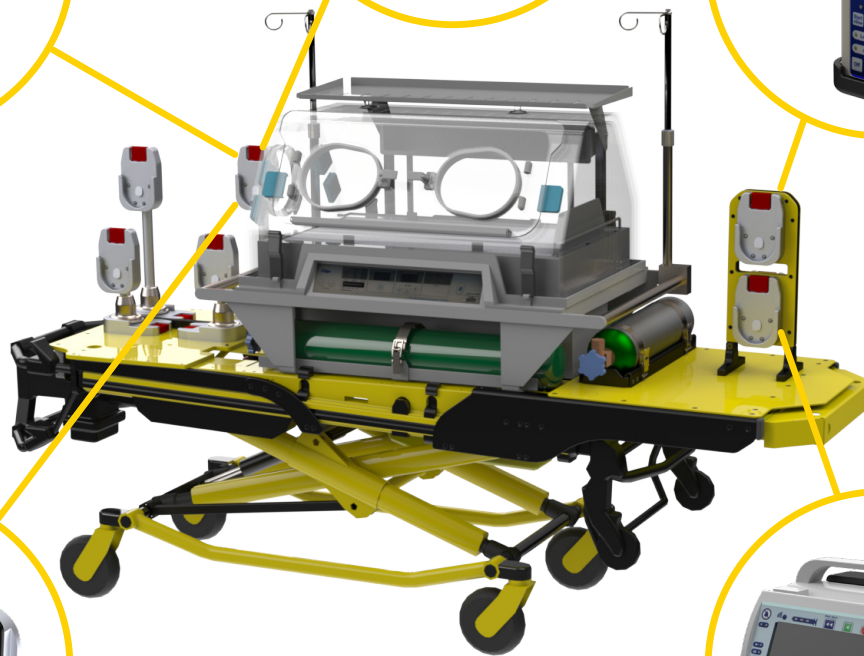

CareFusion Bird®  
Oxygen Blender Bracket



Smiths Medical Medfusion® 4000  
Syringe Pump Bracket



# TECHNICAL SPECIFICATIONS

<b>Product Name</b>	Neonatal Stretcher System
<b>Product Description</b>	Mounting Systems for Medical Devices on the Stryker Power-PRO™ IT Stretcher during Neonatal Critical Care Transport
<b>Compatibility of Stretcher/Equipment</b>	Stryker Power-PRO™ IT with Dräger Isolette® TI500 incubator
<b>Compatibility of Medical Devices*</b>	Customizable Configuration, Several Devices and Accessories Available
<b>Environment of Usage</b>	EMS/CCT (Ground) – Hospital/Clinic
<b>Product Section Material</b>	<ul style="list-style-type: none"><li>• Patient Head Support Structure Side Rails: Aluminum 6061-T6</li><li>• Patient Foot End Support Structure Side Rails: Aluminum 5052-H32</li><li>• Patient Head and Patient Foot End Support Structures: Aluminum 6061-T6</li><li>• Micro Bases: Aluminum 6061-T6</li><li>• Safety MD-Transporter: Aluminum 6061-T6 and Stainless Steel 304</li><li>• Vertical Mounting System: Aluminum 6061-T6</li></ul>
<b>Product Finish</b>	<ul style="list-style-type: none"><li>• Patient Head Support Structure Side Rails: Black-Anodized Finish</li><li>• Patient Foot End Support Structure Side Rails: Black-Anodized Finish</li><li>• Patient Head and Patient Foot End Support Structures: Yellow Powder-Coated Paint</li><li>• Micro Bases: Silver-Anodized Finish</li><li>• Safety MD-Transporter: Silver-Anodized Finish</li><li>• Vertical Mounting System: Black-Anodized Finish and Yellow Powder-Coated Paint</li></ul>
<b>Environmental Regulations</b>	Complying with RoHS 2, RoHS 3
<b>Temperature of Usage</b>	-31° F to 113° F (- 35° C to 45° C)
<b>Certification of Individual Systems**</b>	Designed to comply with: <ul style="list-style-type: none"><li>• SAE J3043 (26 G)</li></ul>
<b>Transportation Compatibility System</b>	<ul style="list-style-type: none"><li>• Safety MD-Transporter (Vertical)</li><li>• Vertical Stretcher Mounting System</li><li>• Micro Bases</li></ul>
<b>Features</b>	<ul style="list-style-type: none"><li>• Allows the Installation of Several Medical Devices</li><li>• Access to all Medical Device Connectors and Accessories</li><li>• No Interference with Incubator Doors for Emergency Operations</li><li>• Allows Rotation of Medical Devices on Safety MD-Transporter</li><li>• Made of High-Density Aluminum</li></ul>
<b>Weight</b>	Between 80 lbs. (36.29 kg) and 85 lbs. (38.56 kg)
<b>System Physical Overall Dimensions</b>	Width 78.0 in, Depth 22.7 in, Height in relation to support structure vary between 22.3 in and 25 in depending on the equipment chosen
<b>Complete System Weight (with brackets, without medical devices)</b>	Contact Technimount EMS for more information on the weight of a specific configuration at <a href="mailto:customerservice@technimount.com">customerservice@technimount.com</a>
<b>Part Number</b>	<ul style="list-style-type: none"><li>• 521-00-SITFP-2 - Patient Head Support Structure</li><li>• 521-00-SITBP-1 - Patient Foot End Support Structure</li><li>• 523-00-SITWAL - Vertical Mounting System</li><li>• 540-00-OL13 - Safety MD-Transporter</li><li>• 544-00-OL08 - Safety MD-Transporter SV (Short Version)</li></ul>

\* Product and company names of medical equipment are trademarks™ or registered trademarks® of their respective holders. Technimount does not have a commercial relationship with these medical equipment companies.

\*\* This certification does not include the Stryker Power-PRO IT stretcher or Dräger Isolette TI500 incubator.

# NEONATAL STRETCHER SYSTEM

## FEATURES

- ✓ Technimount's solution components designed to comply with the highest industry safety standards for transport of medical equipment, SAE J3043
- ✓ Mounting solution allows for the installation of multiple devices and positioning options on the stretcher during transport
- ✓ Customizable solution for EMS/CCT and hospital environments
- ✓ Access to all medical device connectors and accessories without interference
- ✓ High-density aluminum with powder-coated paint and anodized finish

## BENEFITS

- ⊕ Reduces risks of damage to the equipment and ensures the safety of patients and crew during EMS and Critical Care transport
- ⊕ Safe transfer of medical equipment and simplification of inter-hospital transport logistics
- ⊕ Increase flexibility and access to medical devices from all angles complying with protocol in effect
- ⊕ Scalable solution adapted to specific needs
- ⊕ Integrity of all equipment and stretcher respected
- ⊕ Lightweight, robust and resistant to biological fluids and Quaternary-type disinfectants

# YOUR SAFETY IS OUR PRIORITY



## CONTACT US

### TECHNIMOUNT EMS

3770 Jean Marchand Street, Suite 100-C  
Quebec (QC) G2C 1Y6  
Canada  
T: + 1 888-639-2758 (North America)

[Info@technimount.com](mailto:Info@technimount.com)

**SAFETY AND FLEXIBILITY  
WHERE IT MATTERS MOST**

[technimount.com](https://www.technimount.com)



Neonatal Stretcher System, Technimount EMS are Registered Trademarks of Technologies CGC Inc. Product and company names of medical devices are trademarks™ or registered trademarks® of their respective holders.

© 2023 Technimount E.M.S. Holding Inc. - All Rights Reserved - Printed in Canada

TS-NSS-202304EN-02