



**TECHNIMOUNT**  
EMS®

## **BRACKET PRO SERIE® 61 - FL**

MOUNTING SYSTEM FOR AIR TRANSPORT OF THE  
HAMILTON-H900 HUMIDIFIER

---



**SAFETY AND FLEXIBILITY**  
**WHERE IT MATTERS MOST**


# BRACKET PRO SERIE® 61 - FL

## Product Components

The Bracket Pro Serie 61-FL is a mounting solution designed for the Hamilton-H900 humidifier during EMS and critical care air transport.

Equipped with either a bottom micro disc, a back micro disc, or both, this bracket offers quick and easy installation in various mounting configurations, facilitating a seamless continuum of care.

Tested in compliance with the highest industry safety standards, this solution stands as a beacon of reliability and convenience to EMS professionals.

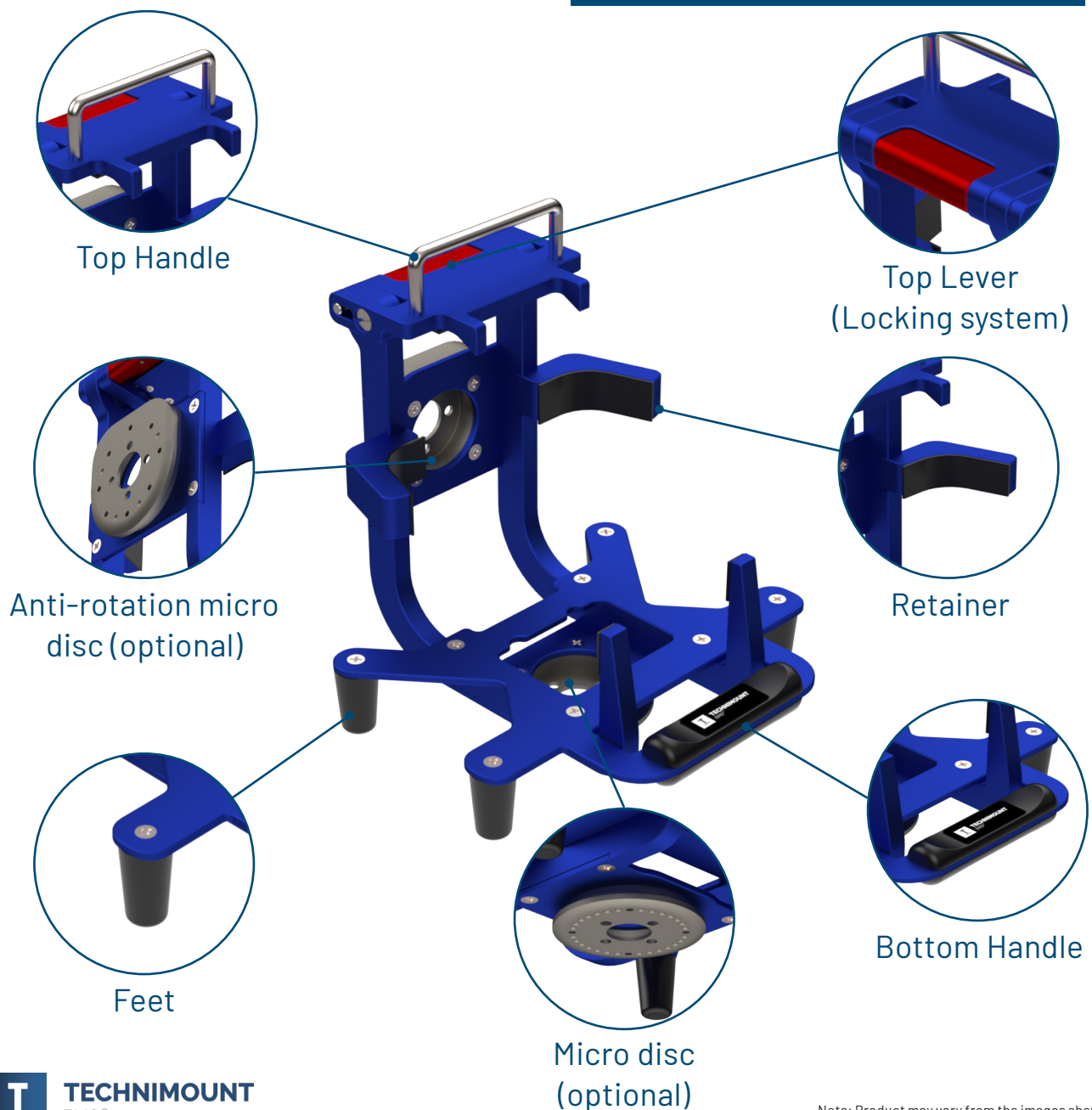

# BRACKET PRO SERIE® 61 - FL

## Highlights

### FEATURES & BENEFITS

- Allows full access to the ventilator's connectors and tubings at all times, helping maintain crucial patient connections.
- Bracket retains and protects the Hamilton-H900 reservoir during transport.
- Equipped with a micro disc (1 or 2) for various mounting configurations on surfaces, stretchers, and wall mounts, accommodating different patient needs during transport.
- Quick release mechanism to easily open the top lever of the bracket to insert or remove the humidifier.
- Durable feet to stabilize the bracket on any flat surface.
- Top and bottom handles allowing efficient installation and removal of the bracket from the base for quick deployment in emergencies.
- Tested in compliance with the highest industry safety standards for air transport of medical devices.

## POSSIBLE CONFIGURATIONS



Micro Base



Safety MD-Transporter

# BRACKET PRO SERIE® 61 - FL

## Technical Specifications

<b>Product Name</b>	Bracket Pro Serie® 61 – FL
<b>Description</b>	Mounting system designed to aid trained EMS and clinical personnel move the HAMILTON-H900 humidifier during EMS and critical care transport
<b>Product Code</b>	Bracket Pro Serie® 61 – FL models: - 710-21-HMH9-BD-FL (back disc) - 710-22-HMH9-2D-FL (back and bottom discs) - 710-20-HMH9-FL (bottom disc)
<b>Operating Environment</b>	EMS/CCT (air and ground)
<b>Compliance</b>	- Tested in compliance with 14 CFR 23.561 & 25.561 & 27.561 & 29.561 (for further details, contact the Regulatory Affairs at customerservice@technimount.com) - Tested in compliance with SAE J3043 and AMD-028
<b>Compatible Stretcher</b>	N/A
<b>Compatible Mounting System</b>	- Micro Base (horizontal) - Micro Base (vertical)
<b>Compatible Medical Devices/ Accessories</b>	HAMILTON-H900 humidifier
<b>Dimensions (W X D X H)</b>	Specifications w/o the medical device: - Back disc: 9.8 in. X 9.6 in. X 11.1 in. (25 cm X 24.4 cm X 28.2 cm) - Back and bottom disc: 9.8 in. X 9.6 in. X 11.1 in. (25 cm X 24.4 cm X 28.2 cm) - Bottom disc: 9.8 in. X 9 in. X 11.1 in. (25 cm X 22.9 cm X 28.2 cm)
<b>Weight</b>	Specifications w/o the medical device: - Back disc: 4.24 lb (1.92 kg) - Back and bottom disc: 4.5 lb (2.04 kg) - Bottom disc: 4.16 lb (1.89 kg)
<b>Composition</b>	- Bracket Pro Serie® 61 – FL: aluminum, stainless steel, plastic, rubber - Micro disc/Anti-rotation micro disc: aluminum
<b>Safe Working Load (SWL)</b>	6.5 lb (2.95 kg)
<b>Operating Temperature</b>	- 31° F to 113° F (- 35° C to 45° C)
<b>Options</b>	N/A